



Other Instrumental Application


Detach the Universal Port from the OCTO Port, in order to insert 15~32mm diameter instruments.

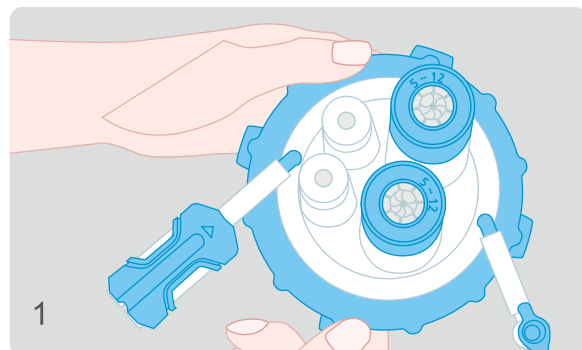


Detach OCTO Port cap from Retractor using Port Cap Lever placed in four ways.

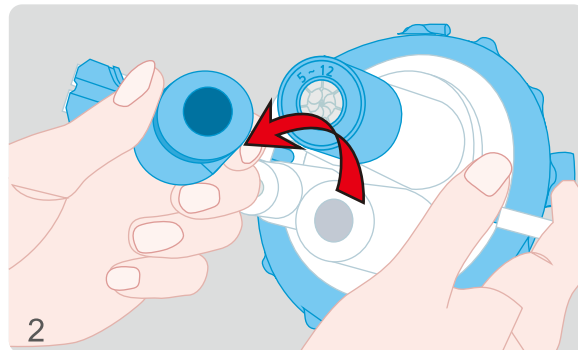

Pull Universal Port strongly twisting it as the bottom of Universal Port is designed to be firmly fixed on silicon cover.

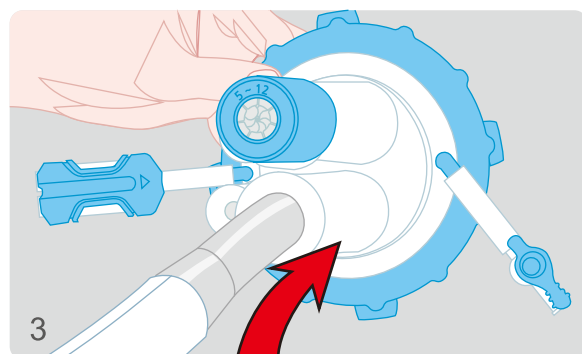

Silicon cover without Universal Port can accept instruments of diameter min.15mm to max. 32mm.

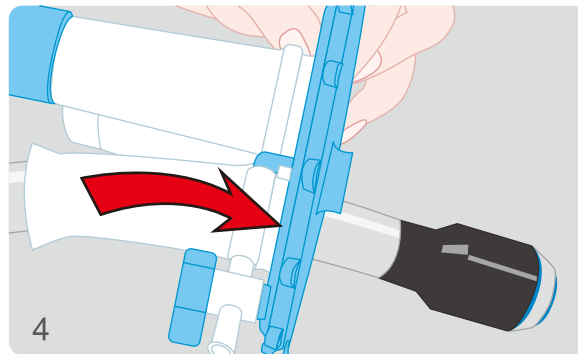

For safer installation, attach OCTO Port cap to Retractor after inserting the instrument correctly.

Universal Port can be installed again on silicon cover in reverse order described above after finishing using the instrument.



The Frankenman Academy is proud to promote several training courses in numerous surgical disciplines with expert faculty to cover both didactic elements and "hands on" skills workshops.

Single port & Reduced port surgery

OCTO™ PORT

OCTO Port allows the user to multitask using its flexible multiple tentacles, responding to diverse situations in surgery.




Single port & Reduced port surgery

OCTO™ PORT



The name OCTO Port originates from Octopus. OCTO Port has many similarities to the Octopus in terms of its appearance and unique features. OCTO Port allows the user to multitask using its flexible multiple tentacles, responding to diverse situations in surgery.



Pure Single Port Surgery

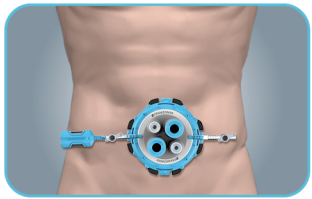

Reduced Port Surgery
Advanced Laparoscopic Surgery

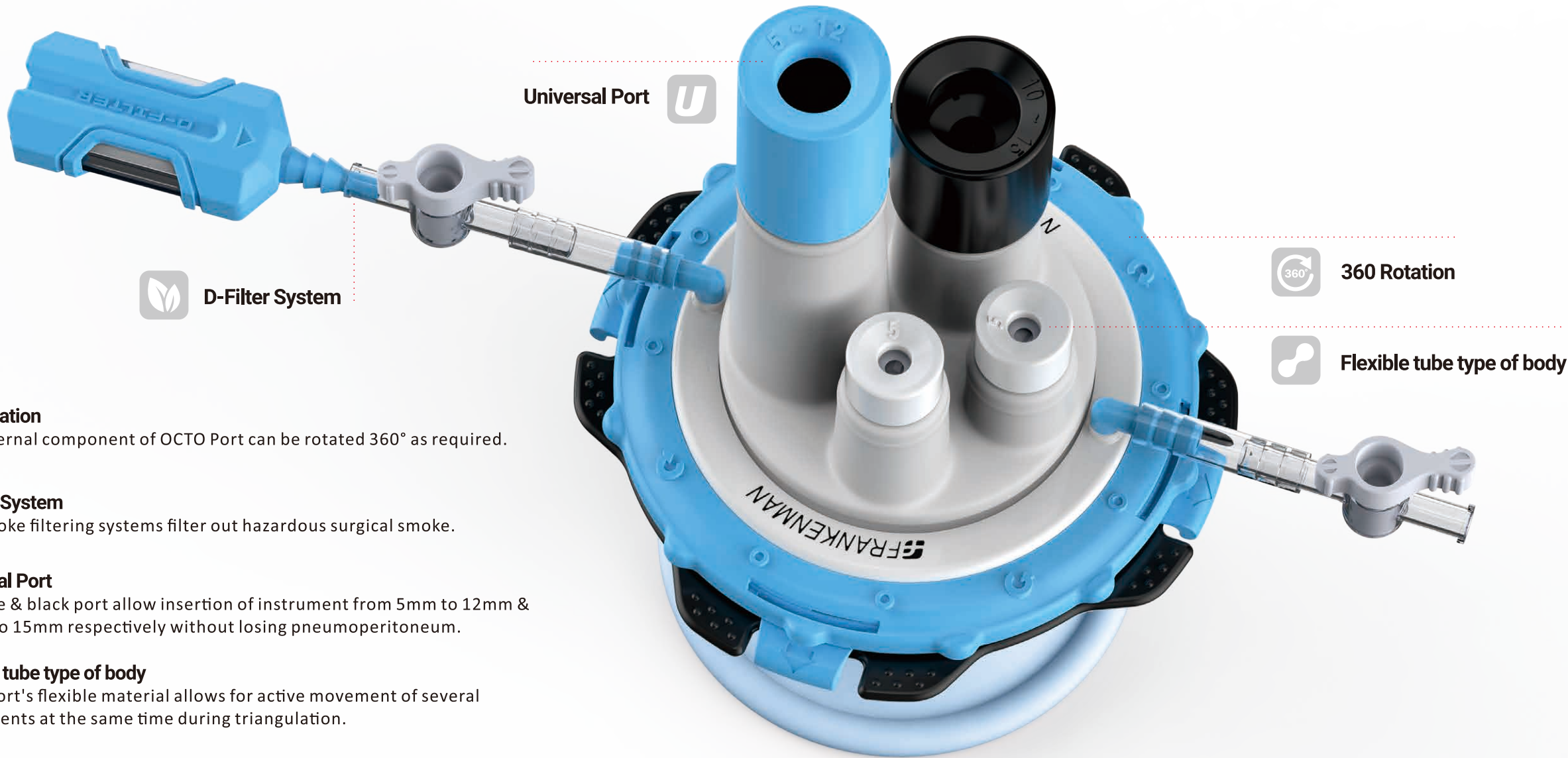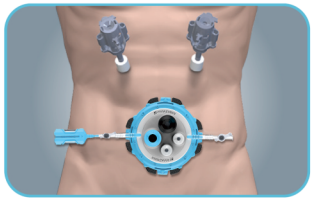

360 Rotation
The external component of OCTO Port can be rotated 360° as required.



D-Filter System
Two smoke filtering systems filter out hazardous surgical smoke.



Universal Port
The blue & black port allow insertion of instrument from 5mm to 12mm & 10mm to 15mm respectively without losing pneumoperitoneum.



Flexible tube type of body
OCTO Port's flexible material allows for active movement of several instruments at the same time during triangulation.

New feature:
Octoport allows product Multitasking using different Valve sizes.

Model Type	Wound Size	5mm Port	5~12mm Port	Quantity
HJN10304A	15~30mm	2	2	6units/box
HJN10504A	35~50mm			

Model Type	Wound Size	5mm Port	5~12mm Port	10~15mm Port	Quantity
HJN10304B	15~30mm	2	1	1	6units/box
HJN10504B	35~50mm				

