


POWERED ENDOSCOPIC STAPLER POWERED BY REUSABLE BATTERY

WITH TISSUE THICKNESS SENSING TECHNOLOGY & SMART FIRING




 **FRANKENMAN**

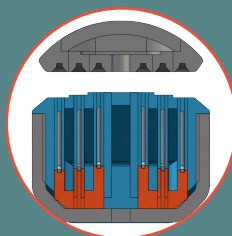

WHEN YOU PUT PATIENTS' HEALTH FIRST,
REMARKABLE THINGS HAPPEN


Frankenman Powered by Reusable Battery. Longer battery life translates directly into a lower environmental impact per use. The more surgical procedures a rechargeable battery can withstand, the fewer batteries will be needed, reducing resource consumption, emissions, and waste production. The reusable battery can be used on average in 10 surgical procedures with approximately 6 applications for each stapler used.

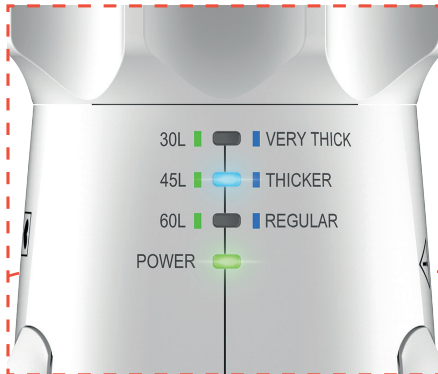


Tissue thickness Sensing Technology: Powered stapler will sense the tissue during the closing and a led will show the different tissue thickness before firing: "Regular", "Thicker" or "Very thick". The integrated chip for Intelligent tissue identification will recognise the right application facilitating the application on that specific tissue.



Smart Firing Speed: Powered stapler will calibrate the speed level to optimise the staple formation, fast speed for the "Regular" tissue, medium speed for "Thicker" tissue, slow speed for "Very thick" tissue.


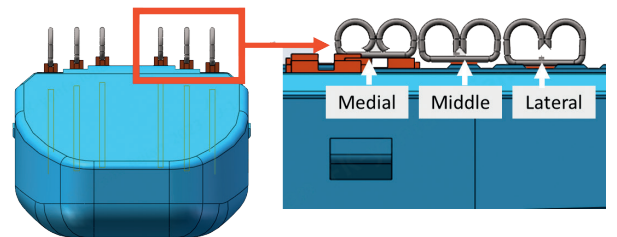


Different height pushers make a different staple forming height, to get a strong anastomosis.

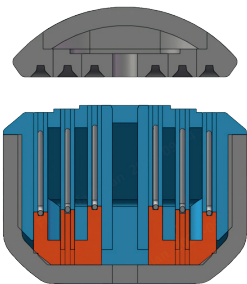


Pushing Staple Technology: Increased balance between hemostasis and blood perfusion in the anastomosis thanks to the new variable compression solution having more margin to maintain closure of the staple leg in all staple lines.

Wavy Grip Solution: The wavy surface of the reload holds the tissue firmly in place, preventing the slippage during the cutting process and retraction of the blade.

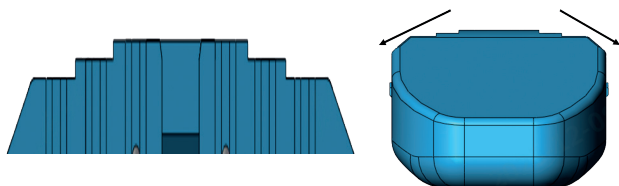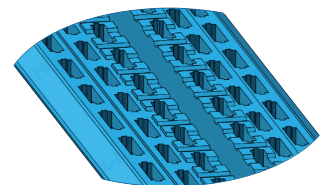

Triple-tiered reload surface: 0,15mm step difference design. The double-sided slanted structure makes the tissue slide on both sides after being compressed for better thickness and proper staple forming.




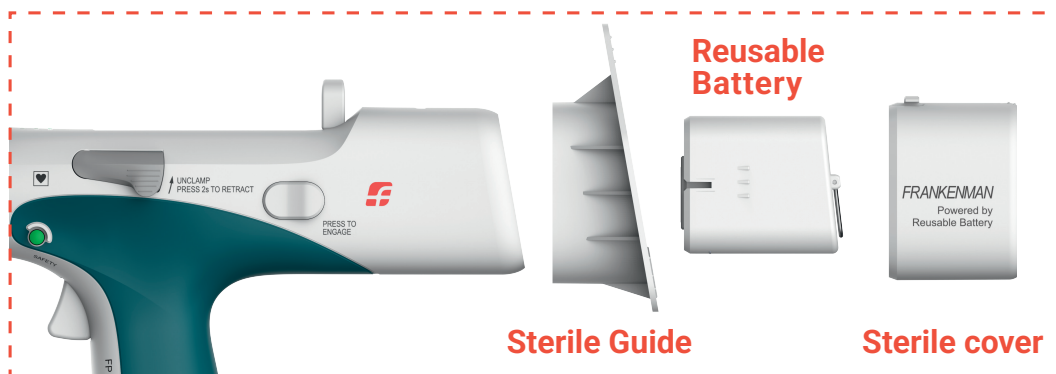

Reusable Battery Insertion:


The reusable battery can be used on average in 10 surgical procedures with approximately 6 applications each with insertion Kit included

Intuitive Control: all controls available on both sides.

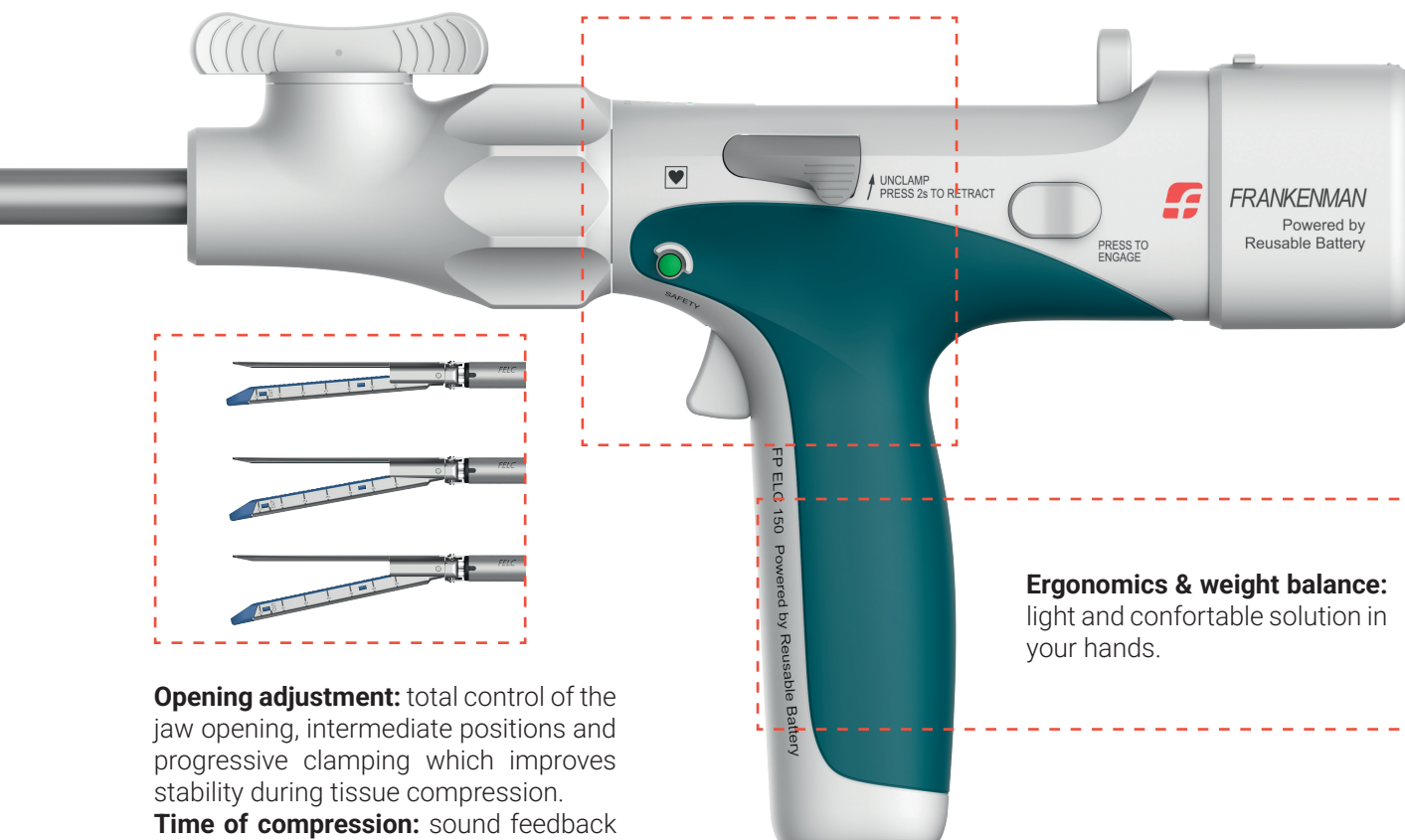


Opening adjustment: total control of the jaw opening, intermediate positions and progressive clamping which improves stability during tissue compression.

Time of compression: sound feedback from 10s to 15s after clamping.

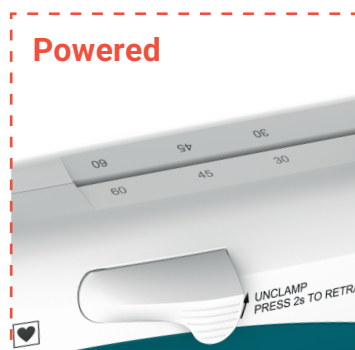

Integrated Overload Protections: an unexpected resistance can be protected automatically. After solving the problem, the stapler can be used for additional firings just replacing the reload. Both powered retraction and manual retraction are allowed.

RB - Reusable Battery: unsterile Reusable Battery pack with charge status.

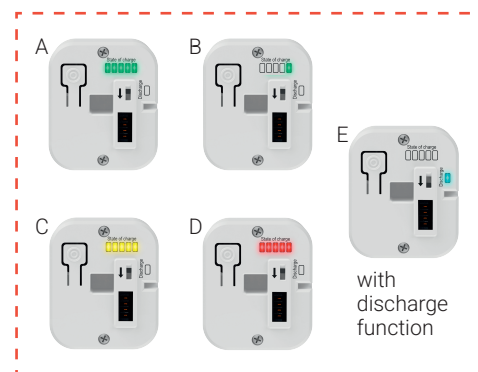

Product Code	Shaft Length (mm)	Working length max (mm)	EXTERNAL DIAMETER (mm)	UNITS PER BOX
FP ELC50	50	270	12.5	3
FP ELC150	150	370	12.5	3
FP ELC250	250	470	12.5	3
RB - Reusable Battery not sterile	—	—	—	3

Model	Colour	Reloading Unit Type	Open Staple Height H/mm	Closed Staple Height H/mm			Degree of Articulation W/°	Staple Quantity	Staple Crown (mm)	Units per Box
NER30-2.0	●	Reversed Stepped Reloading unit	2.0	0.75	0.75	0.75	45	48	3	6
NER30-2.5	○		2.5	1.00	1.00	1.00		48	3	6
NER30-3.5	●		3.5	1.35	1.50	1.65		48	3	6
NER30-4.0	●		4.0	1.60	1.75	1.90		48	3	6
NER30-4.2	●		4.2	1.70	1.85	2.00		48	3	6
NER30-4.8	●		4.8	1.85	2.00	2.15		48	3	6
NER45-2.0	●		2.0	0.75	0.75	0.75		66	3	6
NER45-2.5	○		2.5	1.00	1.00	1.00		66	3	6
NER45-3.5	●		3.5	1.35	1.50	1.65		66	3	6
NER45-4.0	●		4.0	1.60	1.75	1.90		66	3	6
NER45-4.2	●		4.2	1.70	1.85	2.00		66	3	6
NER45-4.8	●		4.8	1.85	2.00	2.15		66	3	6
NER45-5.0	●		5.0	2.15	2.30	2.45		66	3	6
NER60-2.0	●		2.0	0.75	0.75	0.75		90	3	6
NER60-2.5	○		2.5	1.00	1.00	1.00		90	3	6
NER60-3.5	●		3.5	1.35	1.50	1.65		90	3	6
NER60-4.0	●		4.0	1.60	1.75	1.90		90	3	6
NER60-4.2	●		4.2	1.70	1.85	2.00		90	3	6
NER60-4.8	●		4.8	1.85	2.00	2.15		90	3	6
NER60-5.0	●		5.0	2.15	2.30	2.45		90	3	6
NEPR30-2.0	●	Reversed Reloading unit with Curved Tip	2.0	0.75	0.75	0.75	45	48	3	6
NEPR30-2.5	○		2.5	1.00	1.00	1.00		48	3	6
NEPR30-3.5	●		3.5	1.35	1.50	1.65		48	3	6
NEPR30-4.0	●		4.0	1.60	1.75	1.90		48	3	6
NEPR30-4.2	●		4.2	1.70	1.85	2.00		48	3	6
NEPR30-4.8	●		4.8	1.85	2.00	2.15		48	3	6
NEPR45-2.0	●		2.0	0.75	0.75	0.75		66	3	6
NEPR45-2.5	○		2.5	1.00	1.00	1.00		66	3	6
NEPR45-3.5	●		3.5	1.35	1.50	1.65		66	3	6
NEPR45-4.0	●		4.0	1.60	1.75	1.90		66	3	6
NEPR45-4.2	●		4.2	1.70	1.85	2.00		66	3	6
NEPR45-4.8	●		4.8	1.85	2.00	2.15		66	3	6
NEPR45-5.0	●		5.0	2.15	2.30	2.45		66	3	6
NEPR60-4.2	●		4.2	1.70	1.85	2.00		90	3	6
NEPR60-4.8	●		4.8	1.85	2.00	2.15		90	3	6
NEPR60-5.0	●		5.0	2.15	2.30	2.45		90	3	6



The Frankenman Academy is proud to promote several training courses in numerous surgical disciplines with expert faculty to cover both didactic elements and "hands on" skills workshops.

www.frankenman.com | info@frankenman.co.uk | +86-512-68786150

Frankenman International Ltd - 108 Jinfeng South Road, High New District, - Suzhou, 215163, P.R. China